DR-HO'S

neck pain pro

**Less pain and tension
means more activities
and fun in your life!**

Quick Guide | P. 2

Guide rapide | P. 14

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Kurzanleitung | P. 20

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DR-HO'S NECK PAIN PRO

Pain Therapy Device (model P.T.S-IIb)

This device is central to DR-HO'S Neck Pain Pro system. You can use it in combination with the Neck Pain Pro or with the individual sets of body pads. Choose your treatment program and duration easily without having to bend down. Everything is close to hand.

Electrode outlet

Connects to either the Neck Pain Pro or body pads.

Intensity switch

Increase or decrease the intensity of the intended treatment by turning the wheel using your index finger.

Timer button


You can choose from 3 treatment times, namely 10, 20 and 40 minutes. The treatment will end automatically at the end of the chosen program. If you wish, you can extend the treatment time by selecting a program and its duration again, after first switching OFF the device.

Program button

You can choose from 3 programs with different intensity intervals.

Neck Pain Pro

Ergonomic self-contouring design with breathable cover.

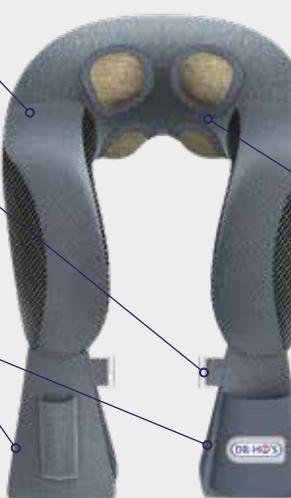Unique built-in electrodes

Clinically proven to help relax tight muscles, improve circulation and relieve pain.

Apply generously DR-HO'S Electro Therapy Conductive Gel over the entire surface of the electrodes before use and reapply every 30 minutes during treatment or if you experience any slight stinging sensations.

Travel clip

Use to secure to luggage or around your neck while travelling.

DR-HO'S Electro Therapy Conductive Gel

This patented gel is developed to moisten the electrodes and the skin, to reduce the amount of contact between the electrode surface and the skin, and to ensure more comfortable and effective treatment.

Mini Spray Bottle

To easily and generously apply the conductive gel over the surface area and edges of the 4 electrodes.

Reusable body pads

The electrode pads, which are available in two sizes, provide effective treatment for specific areas of the body. These zones can be found in the overview on page 6.

IMPORTANT!

DR-HO'S Neck Pain Pro is a medical device class IIA. Please read all guidelines carefully and keep it for future reference.

This complete system (the Neck Pain Pro and the electrode body pads) is intended to relieve pain associated with sore or aching muscles in your shoulders, arms, the back of your neck and other parts of your body, due to strain from exercise or activities of daily living, regular household and work activities.

Neck Pain Pro care instructions

Wipe clean with a moist cloth and allow to air-dry. Do not submerge in water or wash in the washing machine because water may damage the electronic components.

Note:

Do not use the Neck Pain Pro system while sleeping, in the bathroom (when taking a bath or having a shower), driving or operating machinery. Treatment may cause involuntary arm movements and may make you drowsy.

- For optimal results replace the electrode pads after 30 to 50 sessions.
- Electrode pads are intended for personal use and should not be used by others.

Refer to the pad placement chart on how to position the large or small format pads, depending on the area to be treated.

If you wish to order replacement electrode body pads or conductive gel please visit the website of your supplier. If you have any questions, please contact your supplier.

- 1** Insert 2 AAA batteries in the back of the device. It is best to use alkaline batteries.

2 Neck Pain Pro

- A** Plug the wire firmly into the top of the device.

- B** Remove cloths, metal necklaces and jewelry around your neck.

- C** For the best results, apply the conductive gel to your neck and shoulders.

- D** Apply generously conductive gel over the entire surface area and edges of all 4 electrodes.

- E** Rub the gel into the electrodes with your fingers and make sure they are wet.

- F** Place your arms in the arm rests and pull the Neck Pain Pro against your neck. Lean back or lay down.

2 Electrode body pads

- A** Insert the plug from the selected set of electrode pads into the Pain Therapy Device.

- B** To use the pads, remove the transparent protective film from both electrode pads to expose the self-adhesive and sticky gel layer.

- C** Spray a little bit of gel evenly over the surface of the pads, using the handy spray bottle supplied.

- D** Refer to the pad placement chart on page 6. **Note:** The system ONLY works if you attach both pads from the same set to your body.

- E** Attach both pads to the area of the body that you wish to treat. Make sure the pads do not touch each other.

- F** After using the pads, stick the transparent film onto the self-adhesive gel layer again.

- 3** Turn the switch to 'ON'. Then increase the intensity slowly until your muscles start moving of their own accord.

- 4** Press the 'Mode' button to choose from 3 programs with different intensity intervals, A, B or C. A flashing light will indicate the selected mode.

- 5** Press the 'Time' button to select the treatment time, 10, 20 or 40 minutes. Turn the switch to 'OFF' after your 10, 20 or 40 minutes treatment. To continue the treatment first turn the device off and then back on again.

Relax and enjoy!

DR-HO'S NECK PAIN PRO

RELIEVES PAIN & SORE ACHY MUSCLES
at home and at work

INCREASES LOCAL CIRCULATION

SUPPORTS GOOD POSTURE

PAD PLACEMENT CHART

Neck and shoulders

Neck and shoulders

Elbows and forearms

Elbows and forearms

Shoulders

Wrists

Hands


Hips

Lower back

Thighs

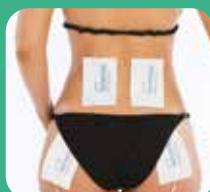Lower back and gluteal muscles

Gluteal muscles and thighs

Knees

Shins

Calves

Ankles

Fast relief for headaches & pain related to sore achy muscles in neck & shoulders.

Disclaimer: DR-HO'S Neck Pain Pro is not designed to diagnose, treat or cure any disease. If you experience any pain or suffer from any health conditions you should consult your doctor for proper diagnosis and treatment.

Precautions: Do NOT use DR-HO'S or any other electronic stimulators if you wear a pacemaker or if you are pregnant, due to unknown safety issues. If you have any questions or concerns, please consult your doctor before use.

“

**Less pain and tension
means more activities
and fun in your life!**

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Guang Zhou Xinbo Electronic Co Ltd.
MedNet EC-REP GmbH